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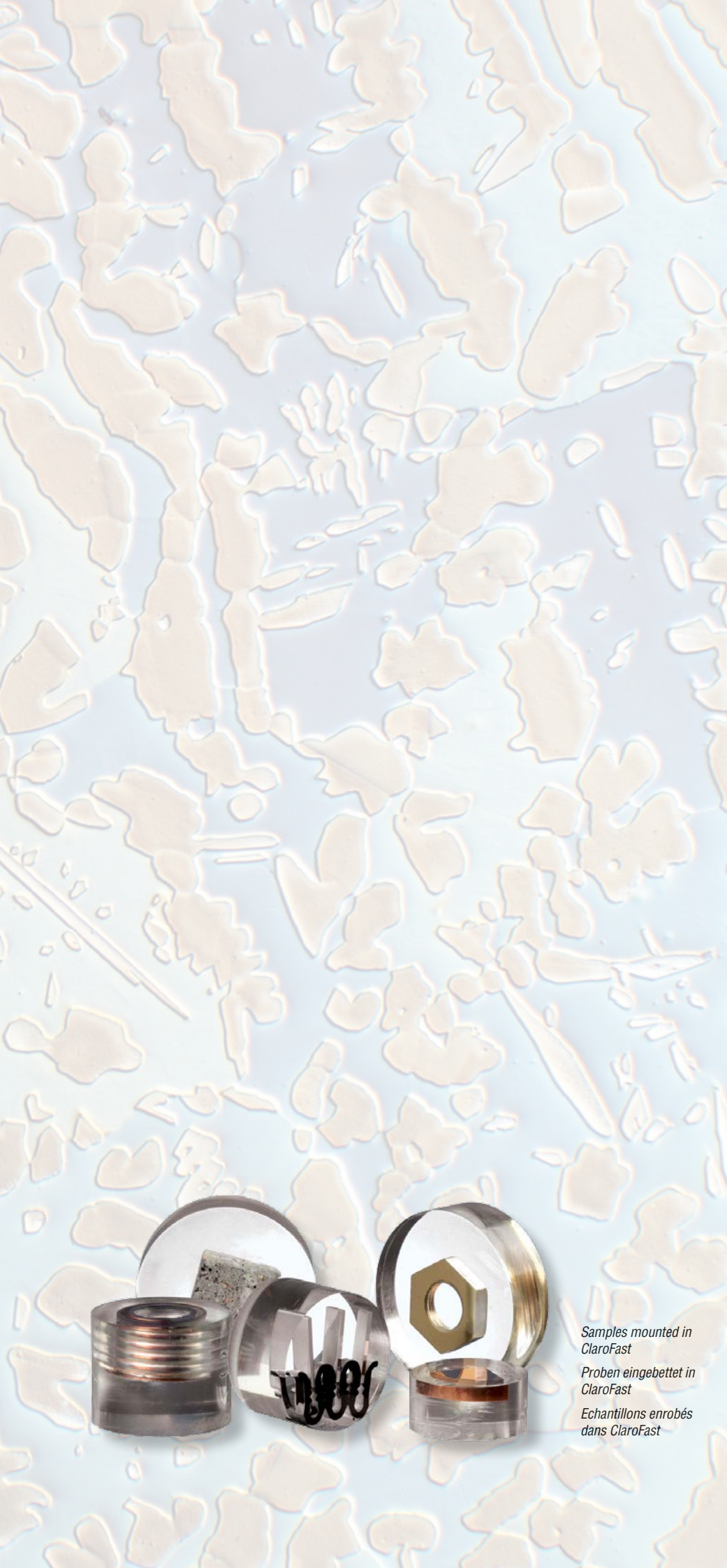
## Table of Contents

	Page	English
Hot Mounting Materials	5	
Mount Release Agent	5	
Accessories	5	
Optimizing the Process Time	6	
Helpful Hints	7	
Mounting Data	8	
Trouble-Shooting	10	

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# Hot Mounting Materials

Product	Application	Specific property	Material	Contents	Cat. no.
ClaroFast	Glass clear mounts Porous samples Surface electrical insulator for ConduFast	Transparent Medium shrinkage	Acrylic Thermoplastic	1 kg 7.5 kg 25 kg	40100055 40100054 40100053
CitoFast	Fast mounting times For soft materials	Very fast mounting. Low process times also when it is used as 'backing' for DuroFast or LevoFast	Acrylic resin with aluminium filler Thermoplastic Grey	1 kg 7.5 kg	40100068 40100069
ConduFast	Electrolytic polishing	Electrically conductive Very low shrinkage	Acrylic resin with iron filler Thermoplastic Black/Grey	1 kg	40100039
DuroFast	Edge retention Planeness	Very low shrinkage Good adhesion Very low removal rate	Epoxy with mineral and glass filler Thermosetting Black	1 kg 7.5 kg	40100044 40100045
LevoFast	Edge-retention for soft to medium hard materials	Very low shrinkage High removal rate	Melamine with mineral and glass filler Thermosetting Light yellow	1 kg 7.5 kg	40100057 40100058
PolyFast	Edge retention	Very low shrinkage Medium removal rate	Phenolic resin with carbon filler Thermosetting Black	1 kg 7.5 kg	40100036 40100037
	Examination in scanning electron microscope	Electrically conductive Low emission in the vacuum chamber			
MultiFast Black	Routine examination Backup resin	Medium shrinkage Medium removal rate	Phenolic resin with wood flour filler Thermosetting Black	2.5 kg 7.5 kg 25 kg 75 kg	40100064 40100065 40100066 40100067
MultiFast Green	Routine examination Colour marking	Medium shrinkage Medium removal rate	Phenolic resin with wood flour filler Thermosetting Green	2.5 kg 7.5 kg	40100078 40100079
MultiFast Red	Routine examination Colour marking	Medium shrinkage Medium removal rate	Phenolic resin with wood flour filler Thermosetting Red	2.5 kg 7.5 kg 25 kg	40100074 40100075 40100076

# Mount Release Agent

Product	Application	Material	Contents	Cat. no.
AntiStick	To prevent the mounts sticking to the rams	Very fine stearate powder	8 g	40300043

# Accessories

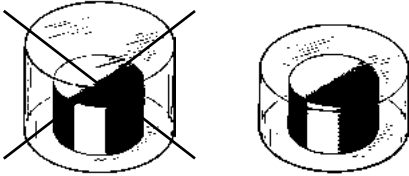
Product	Application	Material	Contents	Cat. no.
Fixation Clips 9 mm	Clip for holding small and thin specimens when mounting	Metal spring	100 pcs.	40300026
Fixation Clips 6 mm		Metal spring	100 pcs.	40300025
EpoDye		Fluorescent dye	20 g	40300002
HandyDoser	Manual dosing unit for hot mounting presses			05736902

# Optimizing the Process Time


When specimens are to be mounted in series it is recommended to optimize the heating and cooling times.

One of the factors limiting both the heating and cooling times is the relatively low heat conductivity of the resins. An efficient way to reduce the process time is to minimize the distance the heat needs to travel through the resin. When mounting metallic pieces (with high heat conductivity), the heating and cooling times may be reduced when:


- The amount of resin is optimized:

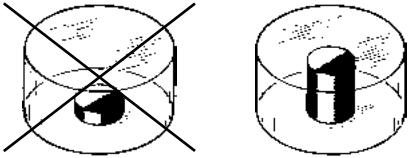
- A relatively high instead of a low sample is mounted:

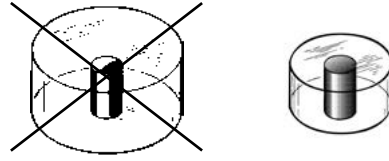- Choose smallest possible mounting cylinder:

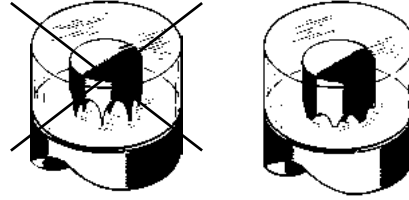- Good thermal contact between metal and ram is ensured:

When optimizing the time it should be reduced stepwise with intermediate inspection of the mount. Insufficient heating and cooling times will result in artifacts (See: "Trouble Shooting").

# Helpful Hints

## ConduFast and ClaroFast

For electrolytical purposes, ClaroFast, although not electrically conductive, can be used in conjunction with ConduFast. By adding a small amount of ClaroFast first, and then completing the required amount with ConduFast, a mount with a conducting body and an insulating preparation surface is formed.

## DuroFast, LevoFast and PolyFast with MultiFast or CitoFast

When using the more expensive PolyFast, DuroFast and LevoFast, a substantial cost saving can be made by adding only a small amount of the desired resin to form the preparation surface, then completing the required amount with the less expensive MultiFast. For fast mounting, use CitoFast as a 'backing' resin. Recommended ratio 1:2 (1/3 DuroFast/LevoFast/PolyFast to 2/3 CitoFast).

## Mount Release Agent

This is recommended to be applied to the mounting rams as a thin layer before the mounting process begins. This prevents the adhesive qualities of the resins from making it difficult to remove the mounts afterwards.

## Distance to Cylinder Wall

The distance between the specimen and the cylinder wall must be a minimum of 3 mm (1/8"), to avoid cracks in the resin. This is especially critical for specimens with sharp corners.

## Small Specimens

Small, thin specimens can be supported during the hot mounting process by the use of Struers Fixation clips. Only metal clips should be used for hot mounting.

## Clean Specimens

To obtain the best results the specimens must be clean, dry and free from grease. If necessary, clean with alcohol or another suitable degreasing fluid.

## Preheating

For porous and/or pressure sensitive specimens, such as minerals, electronic parts etc., it is advantageous to soften the resin by heating, before applying pressure. Preheating is also useful when using thermoplastic resins, particularly ClaroFast.

On CitoPress-15/-30, preheating is available as an automatic programme with the *Sensitive* option. In Sensitive mode, the total heating time is split into two, preheating and heating, where no force is applied during the preheat phase, only during the heating phase.

For metal specimens, we recommend that pressure is applied only for the last minute of the total heating time. If the total heating time is 15 min, preheating should be set to 14 min and 0 bar, and heating to 1 min and 350 bar (ClaroFast).

For mounting of PCB, plastics and other poor conductors, add 1 min to the preheating time (15 min, 0 bar).

## Temperature Sensitive Specimens

The temperature for the mounting process can be reduced to a minimum of 150°C for all resins. This is useful when dealing with temperature sensitive materials. If the temperature is reduced, the recommended heating time should then be increased.

**For very temperature sensitive specimens, hot compression mounting should be avoided. Use Struers cold mounting resins instead.**

## Porous Specimens

Thermoplastic resins, (ClaroFast), penetrate into porous specimens. The best results are obtained by initially preheating the resin (see Preheating).

**For very porous materials, hot compression mounting should be avoided. Use of Struers epoxy cold mounting resins is recommended.**

## Pressure Sensitive Specimens

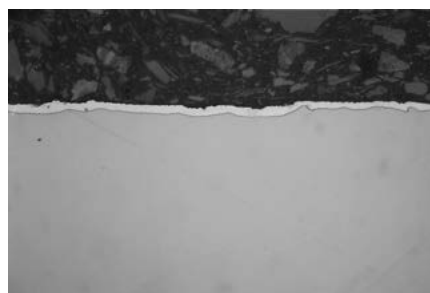
Thermoplastic resins, (ClaroFast), should be used. The best results are obtained by initially preheating the resin (see Preheating).

**For very pressure sensitive specimens, hot compression mounting should be avoided. Use Struers cold mounting resins instead.**

## Specimen Removal

It is crucial that the resin and parameters are correctly matched to each individual specimen. Thermoplastic mounts can be reprocessed.

Specimens mounted in thermosetting resins cannot be reprocessed. They have to be re-mounted. Damage free removal from the completed mount cannot be guaranteed.

*Polished zinc coating,  
Mounted in PolyFast  
1000 x*

# Mounting Data

The recommended mounting parameters appear in the following pages.

## CitoPress-5/-15/-30

The heating and cooling times in the tables refer to the following conditions:

- The heating times in the tables refer to the full process time counted from the start of the process and not from when the pre-set temperature is reached.
- The mounting parameters are based on specimens with a volume approx. 20% of the total volume of the mount. If larger specimens are mounted, less resin should be used.
- If smaller specimens or specimens with low heat conductivity are mounted, the heating and cooling times should be increased. It might also be necessary to increase the pressure to avoid pores in the cured mount.
- A higher pressure is recommended with ClaroFast to avoid "cotton ball" effect.

- The amount of resin in the tables has been adjusted to result in a final height of the mounts of approx. 20 mm (0.8").
- If several resins are combined in one mount, use the process parameters for the resin with the longest times. When using CitoFast as 'backing', use the CitoFast process parameters. For several, or very complicated specimens, add one min. of heating time.

On all CitoPress models, three preset cooling levels are available: High, Medium and Low.

Cylinder dia.	Resin		Heating			Cooling		Time
	Type	Quantity	Time	Temp.	Pressure	Time	Rate	Total time
25 mm / 1"		[ml]	[min]	[°C]	[bar]	[min]		[min]
	ClaroFast	20	4	180	350	6	Low	9.5
	CitoFast	15	2	180	300	1	High	3
	ConduFast	15	3	180	300	1	High	4
	DuroFast	15	3	180	325	2	High	5
	LevoFast	20	3	180	350	1	High	4
	PolyFast	15	3.5	180	325	1.5	High	5
	MultiFast	20	3	180	300	2	High	5

Cylinder dia.	Resin		Heating			Cooling		Time
	Type	Quantity	Time	Temp.	Pressure	Time	Rate	Total time
30 mm		[ml]	[min]	[°C]	[bar]	[min]		[min]
	ClaroFast	20	4	180	350	6.5	Low	10.5
	CitoFast	25	2.5	180	300	1	High	3.5
	ConduFast	20	3.5	180	250	1.5	High	5
	DuroFast	20	3.5	180	325	2	High	5.5
	LevoFast	25	3.5	180	250	2	High	5.5
	PolyFast	20	3.5	180	250	1.5	High	5
	MultiFast	25	3	180	250	2	High	5

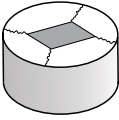
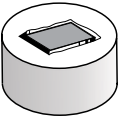
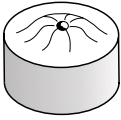
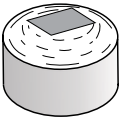
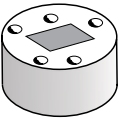
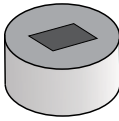
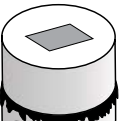
Cylinder dia.	Resin		Heating			Cooling		Time
	Type	Quantity	Time	Temp.	Pressure	Time	Rate	Total time
1¼"		<b>[ml]</b>	<b>[min]</b>	<b>[°C]</b>	<b>[bar]</b>	<b>[min]</b>		<b>[min]</b>
	ClaroFast	20	4	180	350	6.5	Low	10.5
	CitoFast	25	2.5	180	300	1	High	3.5
	ConduFast	20	3.5	180	250	1.5	High	5
	DuroFast	20	4	180	325	2	High	6
	LevoFast	25	3.5	180	250	2	High	5.5
	PolyFast	20	3.5	180	250	1.5	High	5
	MultiFast	25	3.5	180	250	2	High	5.5

Cylinder dia.	Resin		Heating			Cooling		Time
	Type	Quantity	Time	Temp.	Pressure	Time	Rate	Total time
1½"		<b>[ml]</b>	<b>[min]</b>	<b>[°C]</b>	<b>[bar]</b>	<b>[min]</b>		<b>[min]</b>
	ClaroFast	30	4	180	350	7	Low	11
	CitoFast	45	3	180	300	1.5	High	4.5
	ConduFast	35	3.5	180	250	2	High	5.5
	DuroFast	35	4.5	180	350	2.5	High	7
	LevoFast	40	4.5	180	250	2	High	6.5
	PolyFast	30	4.5	180	250	2	High	6.5
	MultiFast	40	4.5	180	250	2.5	High	7







Cylinder dia.	Resin		Heating			Cooling		Time
	Type	Quantity	Time	Temp.	Pressure	Time	Rate	Total time
40 mm		<b>[ml]</b>	<b>[min]</b>	<b>[°C]</b>	<b>[bar]</b>	<b>[min]</b>		<b>[min]</b>
	ClaroFast	35	4	180	350	6.5	Low	10.5
	CitoFast	45	3	180	300	1.5	High	4.5
	ConduFast	40	3.5	180	250	2	High	5.5
	DuroFast	40	4.5	180	350	2.5	High	7
	LevoFast	50	5	180	250	2	High	7
	PolyFast	35	4	180	250	2	High	6
	MultiFast	45	4	180	250	2.5	High	6.5

Cylinder dia.	Resin		Heating			Cooling		Time
	Type	Quantity	Time	Temp.	Pressure	Time	Rate	Total time
50 mm / 2"		<b>[ml]</b>	<b>[min]</b>	<b>[°C]</b>	<b>[bar]</b>	<b>[min]</b>		<b>[min]</b>
	ClaroFast	55	5	180	250	8.5	Low	13.5
	CitoFast	65	3.5	180	250	2	High	5.5
	ConduFast	65	4	180	250	2	High	6
	DuroFast	70	5.5	180	250	3	High	8.5
	LevoFast	75	6	180	250	3	High	9
	PolyFast	55	5.5	180	250	2	High	7.5
	MultiFast	70	6	180	250	4	High	10

# Trouble-Shooting

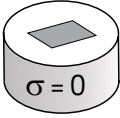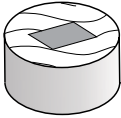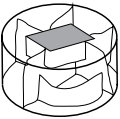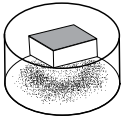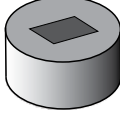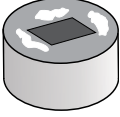General Problems					
<b>Radial Cracking</b>	 <p><b>Cause:</b> Insufficient distance between specimen edge/corner and cylinder wall, or specimen has sharp corners <b>Solution:</b> Increase cylinder diameter or reduce specimen size. The distance between the specimen and the cylinder wall must be a minimum of 3 mm to avoid cracks in the resin. This is especially critical for specimens with sharp corners.</p>	<b>Shrinkage</b>	 <p><b>Cause:</b> Incorrect choice of resin. <b>Solution:</b> Re-mount a new specimen using a resin with a lower linear shrinkage value.</p>	<b>Blistering</b>	 <p><b>Cause:</b> Insufficient heating time. <b>Solution:</b> Increase heating time, or increase process temperature. <b>Cause:</b> Overcured surface. <b>Solution:</b> Decrease process temperature. <b>Cause:</b> Entrapped gas within mount. <b>Solution:</b> Preheat resin.</p>
<b>Bulging</b>	 <p><b>Cause:</b> Insufficient cooling. <b>Solution:</b> Increase cooling time.</p>	<b>Porosity</b>	 <p><b>Cause:</b> Excessive temperature. <b>Solution:</b> Reduce process temperature.</p>	<b>Voiding within large mounts</b>	 <p><b>Cause:</b> Insufficient heating time. <b>Solution:</b> Increase heating time. <b>Cause:</b> Excessive temperature. <b>Solution:</b> Reduce process temperature. <b>Cause:</b> Insufficient force/pressure. Increase mounting force/pressure.</p>
<b>Dull Surface Finish</b>	 <p><b>Cause:</b> Insufficient heating time. <b>Solution:</b> Increase heating time.</p>	<b>Adhesion between mount and rams</b>	 <p><b>Cause:</b> Insufficient application of mould release agent. <b>Solution:</b> Apply mould release agent. This must always be applied to the mounting rams as a thin layer before the mounting process begins. This prevents the resins from sticking to the rams and makes it easier to remove the mounts afterwards. <b>Cause:</b> Insufficient heating time. <b>Solution:</b> Increase heating time. <b>Cause:</b> Excessive force/pressure. <b>Solution:</b> Decrease mounting force/pressure.</p>	<b>Individual Grains Visible on Mount*</b>	 <p><b>Cause:</b> Resin has cured without force/pressure. <b>Solution:</b> Increase force/pressure during heating cycle. <b>Cause:</b> Insufficient heating time. <b>Solution:</b> Increase the heating time and/or temperature.  *Thermosetting resins only</p>

# Trouble-Shooting

ConduFast Problems	
<p><b>Lack of conductivity</b></p> 	<p><b>Cause:</b> No contact with specimen due to use of an excessive amount of ClaroFast, (refer to Helpful Hints).  <b>Solution:</b> Re-mount a new specimen using a smaller amount of ClaroFast than before.  <b>Cause:</b> Insufficient heating time  <b>Solution:</b> Reinsert the mount in the press and reprocess using an increased heating time.</p>
<p><b>Metallic particles in resin removed during the electrolytic process</b></p> 	<p><b>Cause:</b> Insufficient amount of ClaroFast, (refer to Helpful Hints).  <b>Solution:</b> Re-mount a new specimen using a larger amount of ClaroFast.  <b>Cause:</b> Excessive grinding time.  <b>Solution:</b> Re-mount with a new specimen.</p>
ClaroFast Problems	
<p><b>Internal cracking</b></p> 	<p><b>Cause:</b> Excessive cooling rate.  <b>Solution:</b> Reduce the cooling rate</p>
<p><b>“Cottonball” effect in centre of mount</b></p> 	<p><b>Cause:</b> Excessive cooling rate.  <b>Solution:</b> Reduce the cooling rate.  <b>Cause:</b> Insufficient heating time.  <b>Solution:</b> Reduce the physical height of the mount, lower the process temperature and increase the heating time.  <b>Cause:</b> Moist resin.  <b>Solution:</b> Dry the resin by exposing the open container to 30-70°C for 2 hours.</p>
LevoFast Problems	
<p><b>The mounts turn dark after preparation</b></p> 	<p><b>Cause:</b> Insufficient heating time  <b>Solution:</b> Increase heating time and/or temperature</p>
PolyFast and MultiFast Problems	
<p><b>Parts of mount become light or colourless on contact with alcohol</b></p> 	<p><b>Cause:</b> Insufficient heating time.  <b>Solution:</b> Increase heating time and/or temperature</p>

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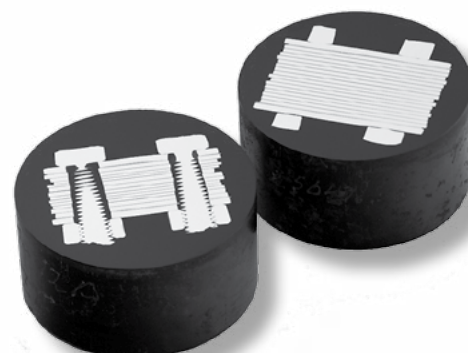
# Struers Application Guide

The following restrictions should be observed, as any violation of the restrictions may cause cancellation of Struers legal obligations:

Struers Application Guide may only be used in connection with Struers products and equipment covered by the Application Guide.

Struers assumes no responsibility for errors in the Application Guide text/illustrations. The information in the Application Guide is subject to changes without notice. The Application Guide may mention consumables, accessories or parts not included in the present version of the Struers product you are using.

The contents of this Application Guide are the property of Struers. Reproduction of any part of the Application Guide without the written permission of Struers is not allowed.

All rights reserved. © 10<sup>th</sup> revised edition, Struers 2011.

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**Struers Ltd.**
**Struers Inc.**