

Product Data Sheet



Struers ApS
Pederstrupvej 84
2750 Ballerup, Denmark

Phone +45 44 600 800
Fax +45 44 600 803
struers@struers.dk
www.struers.com

Product Name **MD/DP-Nap-T**

Cat. no.	40500414	MD-Nap-T, 200 mm (8") dia.
	40500415	MD-Nap-T, 250 mm (10") dia.
	40500416	MD-Nap-T, 300 mm (12") dia.
	40500417	MD-Nap-T, 350 mm (14") dia.
	40500411	DP-Nap-T, 200 mm (8") dia.
	40500412	DP-Nap-T, 250 mm (10") dia.
	40500413	DP-Nap-T, 300 mm (12") dia.

Revision 21-03-2017

Description Short synthetic nap textile with the ability to hold large amounts of suspension and lubricant.

Material Medium-sized regenerated cellulose fibres flocked on woven textile glued on a base textile with an adhesive backing / on a metal plate.

Properties DP-Nap-T to be applied on a preparation disc. DP-Nap-T is self adhesive. MD-Nap-T for magnetic fixation on MD-Disc. Black screen printing on the rear.

Re-usable.

Application Final polishing of all materials. To be used with Struers diamond suspension with grain size 1 – ¼ µm.

Health and Safety

- Whilst the MD/DP-Nap-T is rotating, ensure hands are kept well clear of the edge and out of the splash bowl.
- Do not exceed a rotational speed of 600 rpm.

To the best of Struers' knowledge, the products listed above do not contain any Substances of Very High Concern (SVHC) according to REACH regulations (amended on 20th June 2016) above 0.1%.

According to Article 3(3) of the REACH Regulation (EC) No 1907/2006, this product is considered an article and is defined as *"an object which during production is given a special shape, surface or design which determines its function to a greater degree than its chemical composition"*.

Therefore this product is exempted from CLP obligations and is outside the scope of the UN Globally Harmonized System of Classification and Labelling of Chemicals (GHS).

See Instructions for Use for more information.

Skin protection

Wear protective gloves during polishing.


Respiratory Protection

It is recommended that the polishing machine is adequately ventilated during use of MD-Nap-T / DP-Nap-T

Environmental precautions

Dispose of waste and residues as metal waste in accordance with local authority requirements.

Approved by:

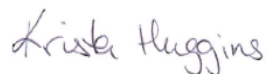

Krista Huggins

Chemical Safety Specialist
Struers ApS.