

Product Data Sheet



Struers ApS
Pederstrupvej 84
2750 Ballerup, Denmark

Phone +45 44 600 800
Fax +45 44 600 803
struers@struers.dk
www.struers.com

Product Name **MD/DP-Nap**

Cat. no.	40500080	MD-Nap, 200 mm (8") dia.
	40500081	MD-Nap, 250 mm (10") dia.
	40500082	MD-Nap, 300 mm (12") dia.
	40500153	MD-Nap, 350 mm (14") dia.
	40500232	DP-Nap, 200 mm (8") dia.
	40500233	DP-Nap, 250 mm (10") dia.
	40500234	DP-Nap, 300 mm (12") dia.

Revision 21-03-2017

Description Short synthetic nap textile with the ability to hold large amounts of suspension and lubricant.

Material Small-sized regenerated cellulose fibres flocked on woven textile with an acrylic adhesive backing / on a metal plate.

Properties DP-Nap to be applied on a preparation disc. DP-Nap is self adhesive. MD-Nap for magnetic fixation on MD-Disc. Black screen printing on the rear.

Re-usable.

Application Final polishing of all materials. To be used with Struers diamond suspension with grain size 1 – ¼ µm and oxide polishing.

Health and Safety

- Whilst the MD/DP-Nap is rotating, ensure hands are kept well clear of the edge and out of the splash bowl.
- Do not exceed a rotational speed of 600 rpm.

To the best of Struers' knowledge, the products listed above do not contain any Substances of Very High Concern (SVHC) according to REACH regulations (amended on 20th June 2016) above 0.1%.


According to Article 3(3) of the REACH Regulation (EC) No 1907/2006, this product is considered an article and is defined as *"an object which during production is given a special shape, surface or design which determines its function to a greater degree than its chemical composition"*.

Therefore this product is exempted from CLP obligations and is outside the scope of the UN Globally Harmonized System of Classification and Labelling of Chemicals (GHS).

See Instructions for Use for more information.

Skin protection	Wear protective gloves during polishing.
Respiratory Protection	It is recommended that the polishing machine is adequately ventilated during use of MD/DP-Nap.
Personal Protection	No requirements.
Environmental precautions	Dispose of waste and residues as metal waste in accordance with local authority requirements.

Approved by:

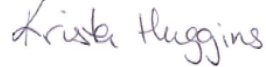

Krista Huggins

Chemical Safety Specialist
Struers ApS.