

# Product Data Sheet



**Struers ApS**  
Pederstrupvej 84  
2750 Ballerup, Denmark

Phone +45 44 600 800  
Fax +45 44 600 803  
struers@struers.dk  
www.struers.com

**Product Name**                      **MD/DP-Pan**

<b>Cat. no.</b>	40500157	MD-Pan, 200 mm (8") dia.
	40500158	MD-Pan, 250 mm (10") dia.
	40500159	MD-Pan, 300 mm (12") dia.
	40500160	MD-Pan, 350 mm (14") dia.
	40500161	DP-Pan, 200 mm (8") dia.
	40500162	DP-Pan, 250 mm (10") dia.
	40500163	DP-Pan, 300 mm (12") dia.

**Revision**                              21-03-2017

---

**Description**                              Polishing cloth. Very high degree of planeness with firm and rigid support.

**Material**                                    Impregnated non-woven technical textile with an adhesive backing / on a metal plate.

**Properties**                                   DP-Pan to be applied on a preparation disc. DP-Pan is self adhesive.  
MD-Pan for magnetic fixation on MD-Disc. Black screen printing on the rear.

Re-usable.

**Application**                              Polishing cloth for fine grinding of soft metals and pre-polishing of hard and brittle materials. To be used with Struers diamond suspension with grain size 15 – 3µm.

**Health and Safety**                      - Whilst the MD/DP-Pan is rotating, ensure hands are kept well clear of the edge and out of the splash bowl.

- Do not exceed a rotational speed of 600 rpm.

To the best of Struers' knowledge, the products listed above do not contain any Substances of Very High Concern (SVHC) according to REACH regulations (amended on 20th June 2016) above 0.1%.


According to Article 3(3) of the REACH Regulation (EC) No 1907/2006, this product is considered an article and is defined as *"an object which during production is given a special shape, surface or design which determines its function to a greater degree than its chemical composition"*.

Therefore this product is exempted from CLP obligations and is outside the scope of the UN Globally Harmonized System of Classification and Labelling of Chemicals (GHS).

See Instructions for Use for more information.

<b>Skin protection</b>	Wear protective gloves during polishing.
<b>Respiratory Protection</b>	It is recommended that the polishing machine is adequately ventilated during use of MD/DP-Pan.
<b>Personal Protection</b>	No requirements.
<b>Environmental precautions</b>	Dispose of waste and residues as metal waste in accordance with local authority requirements.

Approved by:

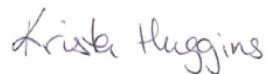

Krista Huggins

Chemical Safety Specialist  
Struers ApS.