

# Product Data Sheet



**Struers ApS**  
Pederstrupvej 84  
2750 Ballerup, Denmark

Phone +45 44 600 800  
Fax +45 44 600 803  
struers@struers.dk  
www.struers.com

**Product Name**                      **MD-Gekko**

**Cat. no.**                      49900047      MD-Gekko, 200 mm (8") dia.  
                                 49900048      MD-Gekko, 250 mm (10") dia.  
                                 49900049      MD-Gekko, 300 mm (12") dia.  
                                 49900050      MD-Gekko, 350 mm (14") dia.

**Revision**                      26-05-2017

---

**Description**                      Adapter for use with SiC Foil or self adhesive backed consumables on the MD-System, for easy removal.

**Material**                      Tinned metal sheet with a high friction surface. A liner protects the surface when the MD-Gekko is not in use.

**Properties**                      The adapter surface has adhesive properties in combination with smooth polymer or metal surfaces. For magnetic fixation on MD-Disc.

Re-usable.

**Application**                      Fixation of non-adhesive polymer or metal-backed as well as self-adhesive backed grinding/polishing consumables.

**Health and Safety**                      -      Whilst the disc is rotating, keep hands well clear of the edge and out of the splash bowl.

To the best of Struers' knowledge, the products listed above do not contain any Substances of Very High Concern (SVHC) according to REACH regulations (amended on 20th June 2016) at or above 0.1%.

According to Article 3(3) of the REACH Regulation (EC) No 1907/2006, this product is considered an article and is defined as *"an object which during production is given a special shape, surface or design which determines its function to a greater degree than its chemical composition"*.

Therefore this product is exempted from CLP obligations and is outside the scope of the UN Globally Harmonized System of Classification and Labelling of Chemicals (GHS).

See Instructions for Use for more information.

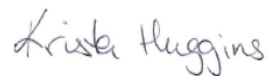
**Skin protection**                      Wear protective gloves during use.

**Respiratory Protection**                      It is recommended that the grinding machine is adequately ventilated during use of MD-Gekko.

**Environmental  
precautions**

Dispose of waste and residues as metal waste or landfill in accordance with local regulations.

Approved by:

A handwritten signature in purple ink that reads "Krista Huggins".

Krista Huggins

Chemical Safety Specialist  
Struers ApS.